



The eSpray is a unique punch face and die wall lubrication system that continuously sprays dry lubricant onto tablet press tooling, delivering a consistent and verifiable amount of lubricant onto the product outer surface.

It is the ideal solution to cope with issues such as capping, sticking, too high ejection force, chemical incompatibility with formulation or poor dissolution.

The system applies only the needed amount of lubricant while vacuuming the excess directly from the nozzle to keep the production environment clean.

It therefore helps to reduce the amount of dry lubricant in the blend, or even eliminate it in the case of new formulations.

The minimal quantity of lubricant sprayed externally onto the tablets can be accurately measured and validated as below the threshold for listing as an ingredient in the formulation.

Recommended in
hygroscopic formulation



Prevent capping and
sticking issues





Key benefits :

- **Accurate spray rate** thanks to high precision feeder
- **Consistent amount** of lubricant sprayed thanks to the venturi pressure control
- Quantity of lubricant sprayed on tablets **below 0,1%**
- High under pressure vacuum integrated to **prevent lubricant recirculation**
- **Retrofittable** to all existing tablet presses
- Available for **single and double** sided tablet press
- **Fast & easy** dismantling and cleaning

Technical data	eSpray
Total Height (A)	1550 mm *
Floor dimension (BxC)	590 mm x 870 mm
Packaging	Wooden Box
Netto weight	180 Kg
Brutto weight	200 Kg
Tablet size	1-28 mm
Spraying range	1,5-15 g/min
Requested compressed air	450 - 600 l/min - Min. 6 bar
Compressed air connection	KS4
Noise level	<80 dB
Power supply	110 V / 230 V \pm 10 %
Frequency	50 Hz / 60 Hz
Non-Metallic Contact Parts	POM C - PU - Silicone
Frame parts	Stainless steel 304
Frame finishing	Brushed 320
Frame mounted on wheels	Yes
IQ-OQ protocol	Available
Material certificates	Available
Hepa Filter	Standard
Available options :	
Dispensing system	For closed container refill

* Depending on options

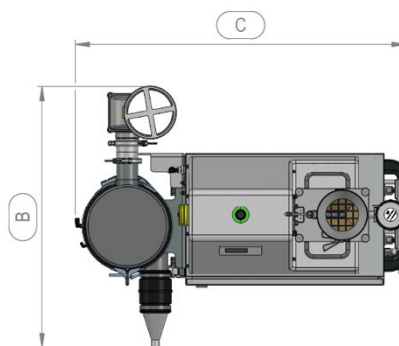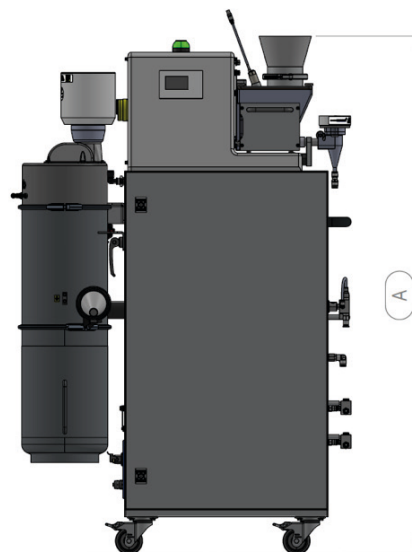

Dust tight dispensing system



Easily retrofittable to all tablet presses



PHARMA  TECHNOLOGY



Rue Graham Bell 8
B-1402 Nivelles, BE
Tel: +32 67 70 13 00
info@pharmatec.be
www.pharmatec.be